OPED bewegt

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OPED gets you moving

We at OPED are experts in medical technology. We develop all our orthopedic rehabilitation products in-house and manufacture them in Germany. Our goal: A mobile patient. **www.oped.de**

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E Open velcro straps for bicep (a) and wrist (b), and then unbuckle belt. (c)
Detach the hook (d) from the arm rest (f).

Open chest buckle (e).
Sit with CAMOshoulder placed on the edge of a table with the straps hanging down.

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E Make sure both shoulders are at the same height.

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E Make sure the arm rest is on the side of your body.

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E Swivel clip (h)
Ensure swivel clip is located over your sternum.

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D Hüftgurt (c)
E Belt (c)
F Sangle de hanche (c)
I Cintura (c)

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E Carefully place your arm in the arm rest (f) and close the velcro straps on your bicep and wrist.

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E Pull the belt (c) around your waist and then close and tighten. The baseplate (i) of the armrest should rest on your hip or pelvic bone.

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E With your good arm, pull the shoulder strap (j) up to rest on the healthy shoulder. The opening on the shoulder pad should face away from the body. Slide the quick release clip (g) toward the belt clip.

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E Attach the hook clip (d) to the arm rest. (f) Open the chest buckle (e) and pull upward until the strap is taut. The swivel slip (h) should always be located over the sternum.

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
E To adjust the arm rest, press the button and set the desired length.

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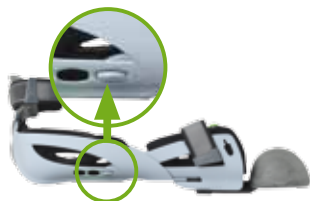
E To adjust the palm rest, squeeze the brackets on the underside of the arm rest and move to the desired position. Palm rest is removable.

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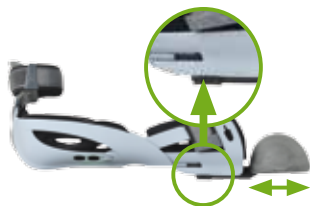E Your CAMOshoulder has already been set to the position that your doctor or physical therapist feels is best for you. Please consult your physician before making any adjustments.

Internal/External Rotation

Internal rotation = rotation of the forearm inwards

External rotation = rotation of the forearm to the outside

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Upon consultation with your
physician/therapist

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E 1. Fixed Position

Levers are next to each other in closed/fixed position.

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E 2. Open ROM

Allow movement by opening levers in 10° increments.

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E 3. Move between fixed position and open ROM

Press the load lever inwards and slide it to the desired position/angle setting. Repeat this process for the other direction.

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Upon consultation with your physician/therapist

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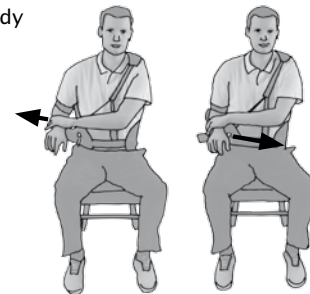
E Open the clip and slide it to the desired angle. To set, press the flap back down. Note: The chest buckle will need to be adjusted to accommodate any changes. (See #5)

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Notice

Cleaning and Maintenance

- Straps and pads: washing machine 40° C (104 F)
- Plastic parts should be hand washed ONLY with a mild detergent
- Air dry only and store in a dry location

Product Features

- Sizes: S/M and L/XL (for each size there is a left and right)

Indications

- Rotator cuff tear
- Tendon ruptures (Supraspinatus, Biceps)
- Endoprosthesis
- AC Joint disorders
- Tuberculum-Majus Fracture
- Dislocations

Safety/Contraindications

There are no known contraindications.

Crashes or impacts of any kind may affect safety. Check your CAMOshoulder immediately for damage. CAMOshoulder should be used exclusively for treatment of shoulder injuries.

Technical data/parameters

Materials

Arm Rest: PP

Small Parts: PA, POM

Straps, Padding: PP, PES, cotton

Screws, Springs: V2A

All parts and components are latex-free.

- Dimensions (without straps):
16.54 in x 8.66 in x 8.27 in"
- Weight:
2.32 lbs

Warranty

The product warranties comply with all legal regulations.

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